



# *SERVICE BULLETIN*

PN: 20812798

**RELEASE DATE:** March 15, 2019

**REVISION DATE(S):** May 13, 2014 (rev A), March 15, 2019 (rev B)

**MODEL:** CAIRE HELiOS Marathon H850


**ISSUE:** Vent Valve Tube Assembly

**NOTES:** The vent valve tube have a threaded connection to the QDV adapter as part of rolling changes to the HELiOS Marathon H850. This is a design improvement to the previous pressed-in connection that will make a more secure connection between the vent valve tube and the adapter. This is not a field replaceable component. This change will not affect the performance of the unit in any way.




The new part number for this assembly is: 20660433

New Design

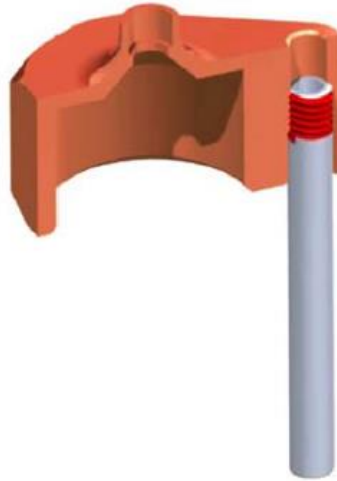

Old Design

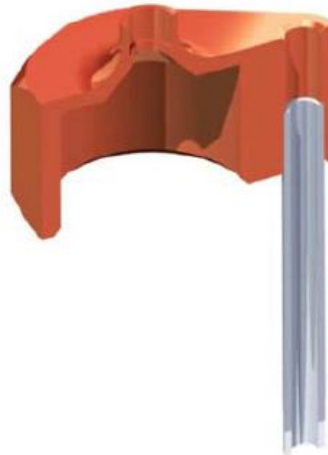

**CONTACT:** For technical questions or concerns, contact Technical Service:

USA: 800-482-2473

[www.cairemedical.com](http://www.cairemedical.com)

Europe: +44 (0) 118 936 7060

[www.cairemedical.com](http://www.cairemedical.com)

For ordering information or general inquiries, contact Customer Service:

USA: 800-482-2473

[www.cairemedical.com](http://www.cairemedical.com)

Europe: +44 (0) 118 936 7060

[www.cairemedical.com](http://www.cairemedical.com)

*Copyright © 2019 CAIRE Inc. CAIRE Inc. reserves the right to discontinue its products, or change the prices, materials, equipment, quality, descriptions, specifications and/or processes to its products at any time without prior notice and with no further obligation or consequence. All rights not expressly stated herein are reserved by us, as applicable.*